

Certification



State Food and Drug Administration



NRTL



CE (MDD)



ISO13485 / ISO9001

Sleeves

Picture	Sleeve Type	Order NO.	Chamber
	Arm Sleeve	220041	4
		220081	8
	Leg Sleeve	220042	4
		220082	8
	Waist Sleeve	220083	4
		220083L	
	Full Body Sleeve	220016FB	8
	Calf Sleeve (DVT Sleeve)	220052	2
	Foot Sleeve (DVT Sleeve)	220053	1



Guangzhou Longest Medical Technology Co., Ltd.

Address: 301 & 401 of Building 2 & Building 3, No.96, Chuangqiang Road, Ningxi Street, Zengcheng District, Guangzhou, Guangdong Province, 511399, China

Tel: 86-020-66353999 **Fax:** 86-020-66353920

E-mail: export@longest.cn **Website:** www.longestmedical.com



LGT-2200L

Compression Therapy Device



Health Life Longest Care


Longest



LGT-2200L is sequential pneumatic pressure system that through program controlled inflation and deflation of the sleeves covered in arms, legs, waist, and back to realize therapeutic effect


⊕ Mechanism and functions

The inflation and deflation of the sleeves enhances the backflow of venous blood and lymph, which can improve the blood circulation of diseased areas. This will lead a series of beneficial effect such as improving the function of peripheral blood vessels, removing the edema, accelerating the healing process, and preventing Muscular Dystrophy.

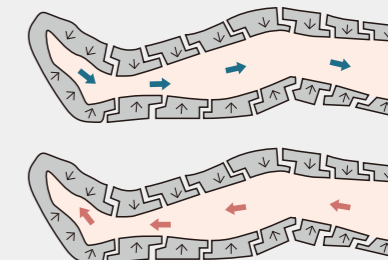


Features

- ⊗ Overlapping air chambers generate a gliding pressure wave with high efficiency
- ⊗ User-friendly interface, just a few steps to finish your setting
- ⊗ Gentle inflation and deflation
- ⊗ Real-time pressure, time, and inflation part monitoring (visible in the interface)
- ⊗ Durable & easy-use sleeves
- ⊗ Pressure unit switchable between mmHg and kPa

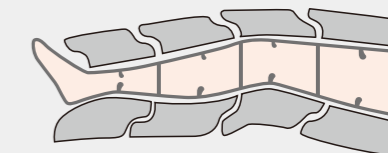

Overlapping Air Chamber

Four chambers modes

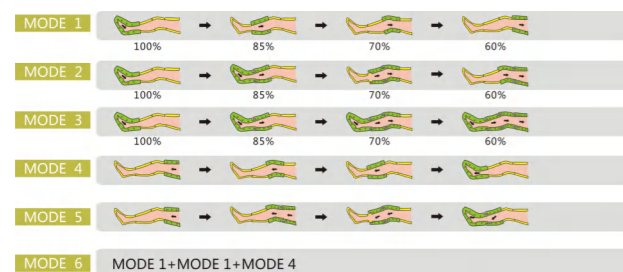

⊕ Specifications

Channel	2
Sleeve chambers	1 / 2 / 4 / 8
Compression modes	18 (Including 3 DVT modes + 6 four chambers modes+ 9 eight chambers modes)
Time range	1 - 99 mins
Pressure range	40 - 190 mmHg
Interface	5.7" touch screen
Dimension	393 (L) x 366 (W) x 194 (H) mm
Weight	10kg

⊕ Application

- ⊗ Prophylaxis of Deep vein thrombosis (DVT) (Post-stroke deep vein thrombosis, Post-traumatic deep vein thrombosis)
- ⊗ Prophylaxis and treatment of edema (Chronic venous edema, Lymphedema, Post-mastectomy lymphedema)

